




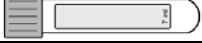


ORDER OF BLOOD DRAW

TUBE COLOR	COLLECTION TUBE	MIX BY INVERTING IMMEDIATELY
Blood Culture Bottles O2/ANO2	Blood cultures	8 – 10 times
LIGHT BLUE 	Citrate	3 – 4 times
GOLD/SPECKLED 	Serum separator tube	5 times
DARK GREEN  LIGHT GREEN 	Heparin Tube Plasma separator tube	8 –10 times
LAVENDER 	EDTA	8 – 10 times
GRAY 	Fluoride tube	8 – 10 times